

**SB<sup>®</sup> 755**



When comfort comes first



Combining comfort, ease of handling and safety, the *Invacare SB 755* is ideal for domestic and community care environments.

The *SB 755* offers optimal comfort with an ideal seating positioning to minimise shear and friction. Carers can adjust the bed to meet individual user requirements.

Thanks to a unique integrated feature, the *SB 755* can be adapted to taller users as both the head and foot of the bed can be extended.

All these features provide real benefits to users and carers in terms of comfort and general well-being.



### Your comfort is our focus

The *SB 755* comes with several features to ensure an excellent seating position and to reduce shear and friction during movement. These include improved auto-regression, optimised length of mattress platform sections and a leg section that can be lowered below horizontal.

Invacare strives to cater for all patients and so offers a length extension of 100 mm at the head and foot of this model to meet the needs of taller users.



Now available with a Light Oak colour option\*

### Ease of handling

No part of the *SB 755* weighs more than 25 kg and dismantlable scissors are available as an option.

Mounting and dismantling the bed has been made particularly easy to allow for greater ease of transport and storage. The bed can be stored on its base and rolled directly from van to bedroom. Assembly and transportation can be carried out by one person quickly and with minimum effort.


### Safety is key

As always at Invacare, safety is our key concern. Side rails in both the standard bed length and the extendable version meet the increased strength demands of the IEC 60601-2-52 standard for beds, without compromising the ergonomics or working conditions of the carers.

### Features and Options



Auto regression – improved functionality. Increased comfort when activating backrest, minimising shear/friction forces.

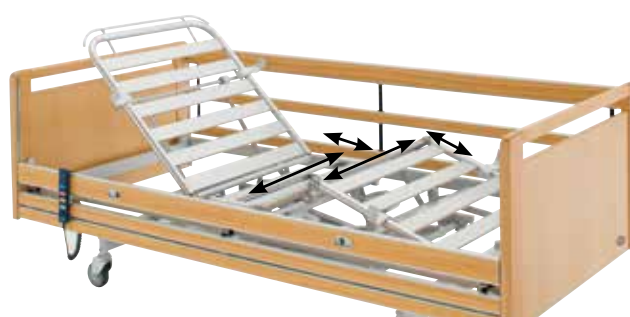

Optimised dimensions of mattress platform according to anthropometrical measurement, for increased comfort.



Extendable mattress platform. Additional 100 mm at head and foot. Bed easily extends with 50, 100, 150 or 200 mm to suit user needs.

---

## Features and Options

---



### Line Side rail

Full-length side rail meets the strength and safety demand of the IEC 60601-2-52 bed standard. The Line side rail is lightweight and easy to use.



### Line Extendable side rail

Full-length side rail, which can be adjusted incrementally from 2000-2200 mm.



### Removable panel

On the Victoria bed end, a removable panel allows easy access from the head and foot of the bed even with the full-length side rails mounted.



## SB 755 Detachable




Bed ends can be dismantled without tools.

**Need to move the bed from one place to another?**  
**No lift?**  
**Narrow staircase?**  
**...No problem!**

The **SB 755 Detachable** is designed to overcome all these challenges.

This bed is so easy to mount and dismantle that it can be carried out by one person; our unique dismantlable system is designed for minimum effort - no need to worry about losing screws or parts!

Designed with parts that weigh no more than 25kg each, the **SB 755 Detachable** is a novel design that is unique in the market. Such easy handling ensures that the bed is easy to carry - to whatever floor.

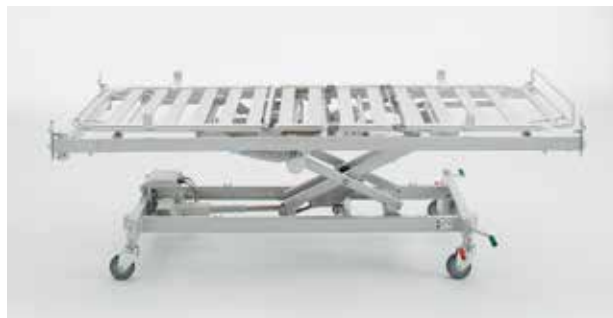

The seat section opens and the two mattress support parts are divided. The mattress platform is detached from the base with a quick release mechanism.



Mattress platform is effectively released and unlocked in the seat section.



The scissor is released from the base with the quick release handle.



Activate the quick release and the scissor is detached. The handle makes it very easy to dismantle.



Bed parts, mattress and accessories are stored directly on the base with special hooks; storage and transport are facilitated and the bed takes up less space.

## Features and Options



**Victoria bed end**

Detachable panel - Available in beech. Compatible with full length side rails.



**Susanne bed end**

Available in beech. Compatible with full length side rails.



**Anita bed end**

Available in beech. Compatible with full length side rails.



**Cover for wooden side rail**

Padded cover for wooden side rails. Cover remains in place when side rail is lowered.



**\*Olivia bed ends**

Can be used in combination with wooden side panel. Available in light oak only.



**\*Wooden side panel**

Can be used in combination with Scala Basic Plus 2 side rails. Available in light oak only.



**Several hand controls:**

*Invacare Soft Control* or standard with "green" button function.



**Lifting pole**

The handle can be adjusted both in height and depth.




**Swivel support handle**

The swivel support handle provides good support for the client getting to a standing position.




**Mattress support**

Mattress support with aluminium slats.

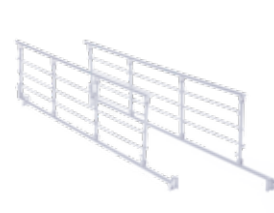


**Invacare Scala Basic 2**

Foldable steel side rail.\*

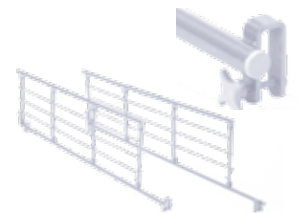


**Invacare Scala Basic Plus 2**

Foldable steel side rail with plastic inserts at the mounting brackets.\*

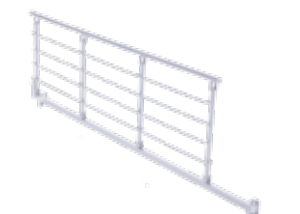

**Invacare Scala Medium 2**

Foldable steel side rail for higher mattresses.\*

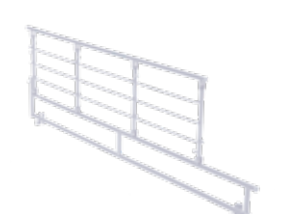

**Invacare Scala Decubi 2**

Foldable steel side rail for extra high mattresses.\*



**Central brake**

Available as standard.



**Lowering kit**

Using castors of 75 mm creates a lower height of 250 mm.

\*Compliant with the new standard for beds IEC 60601-2-52.



**SB755 Wide -Care solutions for plus sized people**

Comprising all the features of the **SB755**, the **SB755 Wide** is available in two sizes with a mattress platform width of 105cm and 120cm. The **SB755 Wide** supports plus sized clients who require additional space for comfort and positioning. The 105cm version is particularly beneficial for carers as it reduces the risk of back injury caused by over stretching across the bed. The **SB755 Wide** has a weight limit of 200kg.

For more information, please see **SB755 Wide** brochure

**Technical data**

For more comprehensive pre-sales information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk)

	Width outside Width inside	Length outside Length inside	Height adjustment	Mattress support dimensions	Under bed clearance	Backrest angle	Thigh angle
<b>SB 755</b>	950/1000 mm 850/900 mm	2130 mm 2080 mm Mattress: 2000mm	350 - 820 mm Total height: 1600 - 2070 mm	800 - 260 - 260 - 680 mm	170 mm	0 - 73°	0 - 28°
<b>SB 755 Wide 1050 mm</b>	1150 mm 1050 mm	2130 mm 2080 mm Mattress: 2000mm	350 - 820 mm Total height: 1600 - 2070 mm	800 - 260 - 260 - 680 mm	150 mm	0 - 73°	0 - 28°
<b>SB 755 Wide 1200 mm</b>	1300 mm 1200 mm	2130 mm 2080 mm Mattress: 2000mm	350 - 820 mm Total height: 1600 - 2070 mm	800 - 260 - 260 - 680 mm	170 mm	0 - 73°	0 - 28°

	Minimum angle below horizontal	Legrest lift	Angle below/ above frame	Max. user weight	Heaviest part of the product	Washable in washing machine
<b>SB 755, SB 755 Wide</b>	0 - 16°	0 - 18°	below: 25° above: 14°	SWL 235 kg (37 stone) 200 kg (31.5 stone)	25 kg	



Bariatric use



**NB**  
The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

**3 year warranty**  
Complies with IEC 60601-2-52

**Frame colour and wood work**

Aluminium grey

Beech

Light oak



The *Invacare SB755* is part of the Safe Patient Handling range.



safepatienthandling

Invacare Limited  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail: [ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

